


Micron Marker

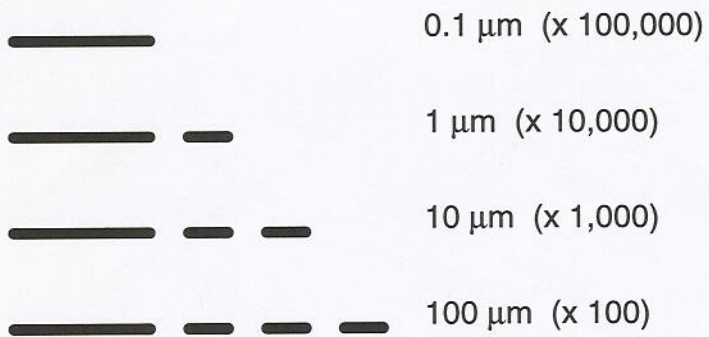

Photography

- 1) Check camera focus. Film ASA 125, use aperture f4.
- 2) Set SCAN MODE to RAPID, Reduced Image area
- 3) Turn CONTRAST and BRIGHTNESS fully anticlockwise.
- 4) Slowly increase BRIGHTNESS to give reading of 10 - 15 on Photo Meter
- 5) Slowly Increase CONTRAST to increase reading further to about 45
- 6) Set SCAN MODE to PHOTO
- 7) Open camera shutter, and clamp it open
- 8) Start Photo Scan
- 9) Release camera shutter, wind on film.